


HIM main chamber vent

Wednesday, September 8, 2021 2:34 PM

1. In case of sample holder not secured




- Need to vent the main chamber and retrieve the holder from the other side
- a. Keep the transfer door open

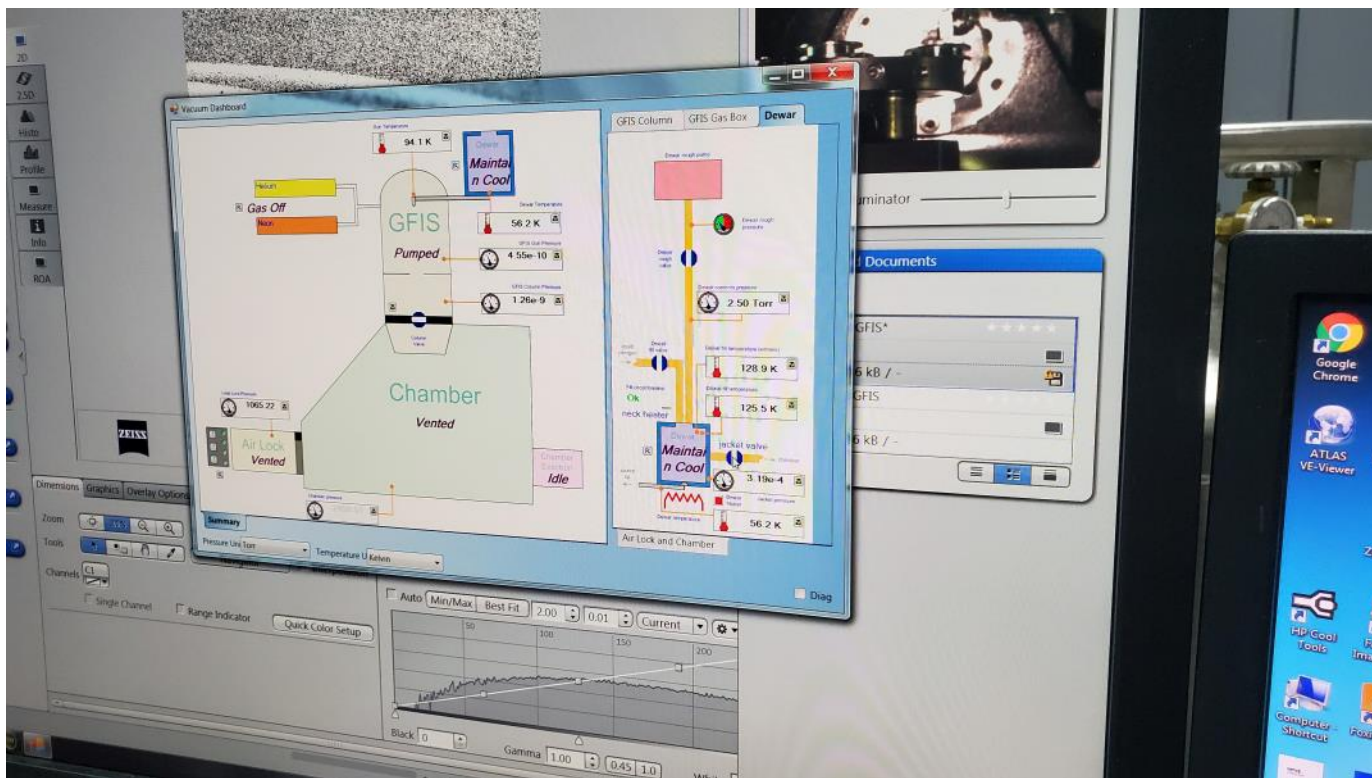

b.

- c. In the vacuum status software, under Dewar tab, make sure the jacket valve is closed
- d. In the vacuum status software, right click Main chamber -> vent
- e. Wait till both main chamber and loading dock vent to 1000+ Torr
- f. In GFIS Gas Box tab, turn off the air stable desk

Open main chamber from the other side

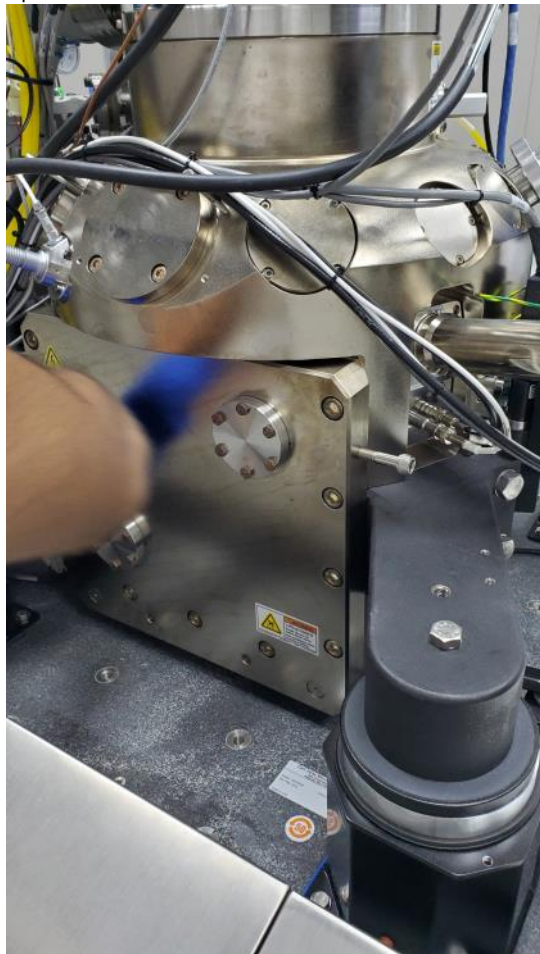

Unscrew the screws, make sure no cables in the way of the chamber door. Place screws on

clean tissue

Open the chamber door, carefully retrieve the sample holder. Make sure not touching ET detector filament or source tip

Close chamber door. First screw on screws loosely by hand, then tighten a few rounds with screw driver. Do not screw too tight.

Load the sample holder from the loading chamber. Because all chambers are vented, you need to hold the loading door when putting the holder back into main chamber

Pump the main chamber from software.

Loading door will close automatically

Turn on table stabilization

Wait till main chamber pump to below 10^{-2} Torr

Pump overnight

Check vacuum the next day and check ion source

Delete user data

D://user data/manohar/HIM/to_drive/service alert.pdf

Follow the instructions to delete user data

In service machine -> hardware communication

Anything wrong can restart zen/server/computer

Power outlet try off and on GFIS high voltage power supplies